



# ULTRAHEALTH™

Experience the complete solution for irregular corneas  
with the **most advanced hybrid lens**



**Long-Term  
Corneal Health**

**Superior  
Vision**


**All-Day  
Hassle-Free  
Wear**

**Intuitive,  
Dependable,  
Efficient Fitting**

# The UltraHealth Advantage

UltraHealth is the most advanced hybrid contact lens with...

## Aspheric Reverse Geometry Vaulted Lens Design

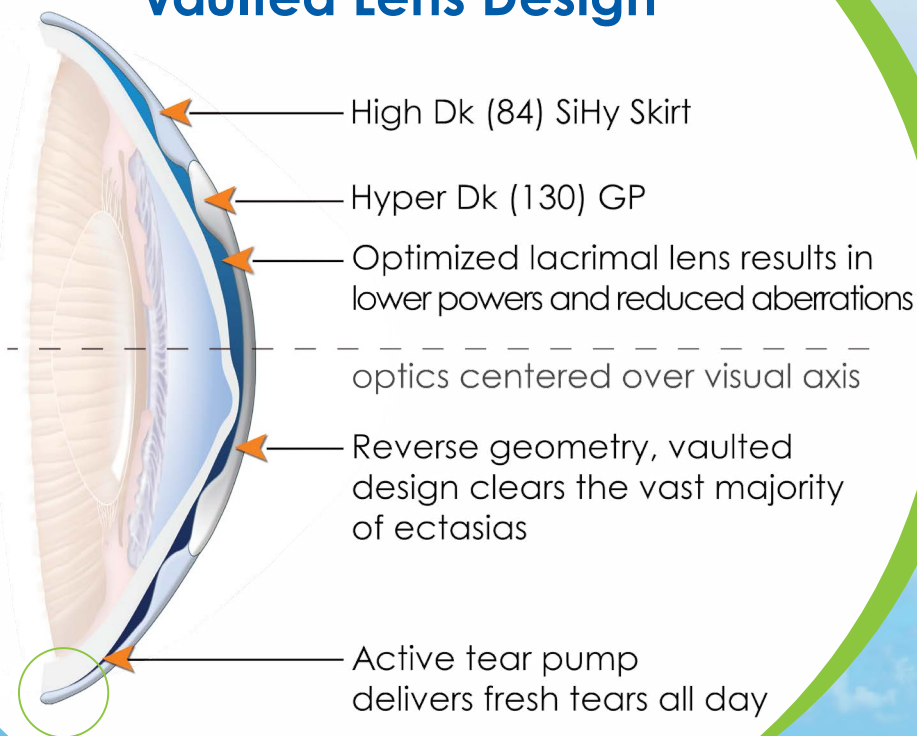

### Long-Term Corneal Health

All-day tear exchange and ultra-high Dk materials

### Superior Vision

GP vision combined with the miracle of centration


### All-Day, Hassle-Free Wear

No more hassle from lenses popping out or from mid-day cleaning

### Intuitive, Dependable, Efficient

Fitting is easy to learn delivering consistent, dependable outcomes.

UltraHealth - OCT Image

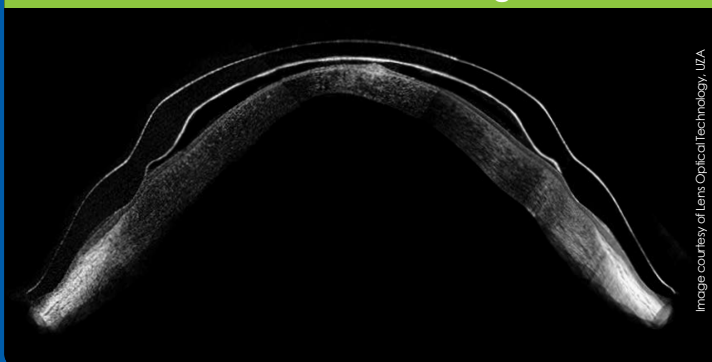

## UltraHealth fitting is:

- Easy to Learn
- Straightforward & Intuitive
- Consistent & Dependable

"Patients can wear their lenses comfortably for 14-16 hours."

— D. Roddy, O.D.

# ge for irregular corneas...

## Long-Term Corneal Health

- Active Tear Pump
- High Dk soft skirt
- Hyper Dk GP

## Superior Vision from 2nd Generation Design

- Vaulting the cornea with centered GP optics provides more reliable and consistent vision
- First hybrid with **new GP design**: reverse geometry, vault system with twice the junction lift
- First irregular cornea hybrid with **new skirt design**: a new outer landing zone for increased comfort
- Enhanced tear exchange from new design

## Hassle-Free

- Soft lens-size makes them more familiar to handle
- Minimal dirt and debris entrapment
- Lenses won't "pop out"
- Wear all day
- No need to remove mid-day for lens build up



## Who is a good candidate for UltraHealth?

Patients with irregular corneas looking to enhance the quality of every day life. This includes young and old, new and existing wearers, GP intolerant patients and soft lens wearers in need of clearer vision. Consider UltraHealth for any patient seeking optimal vision, comfort and corneal health.

## Patient Candidates

- Irregular Corneas
- Emerging to advanced Keratoconus
- Ectatic Corneas
- Intacs®
- Post Corneal Crosslinking
- Post Pk
- Post Rk
- Post Lasik

# The UltraHealth Advantage

## UltraHealth for irregular corneas...with **healthy, hassle-free vision**

Help your patients eliminate the hassles of GP and scleral lenses. Avoid the inconvenience of GP lenses that pop out unexpectedly or debris that lodges uncomfortably under the lens. No need to compromise on comfort. Eliminate patient worries about insertion and removal of large scleral lenses, as well as the hassle of needing to remove, clean and reinsert lenses mid-day.

Eye Health	UltraHealth 15 out of 15 stars	GP's 8 out of 15 stars	Sclerals 7 out of 15 stars
High Dk Materials	★	★	★
Active tear pump under lens	★	★	
Comfort more similar to a soft lens	★		★
Lens vaults the corneal irregularity	★		★
Vision Comparison			
Gas permeable Optics	★	★	★
Optics center consistently over visual axis	★		★
Tear lens less than 100μ	★	★	
Hassle Free			
Lens same or smaller size than soft lens = easier handling	★	★	
Lens removal without special tools	★	★	
Lens does not pop out unexpectedly	★		★
Debris does not freely lodge under the lens	★		★
No need for mid-day lens removal due to build-up	★	★	
Intuitive, Dependable, Efficient			
Majority of patients fit in 3-4 visits	★	★	
Fitting process is easy to learn	★		
Fitting is intuitive and delivers consistent and dependable outcomes.	★		

### Lens Materials & Parameters

#### Lens Materials

- **130 DK** rigid gas permeable center
- **84 DK** silicone hydrogel soft skirt
- Class II UV blocker: >80% of UVA and >95% of UVB radiation

#### Diameter

14.5mm

#### Skirt Curves

8.7 Flat 2, 8.4 Flat, 8.1 Medium, 7.9 Steep

#### Lens Powers

##### +10.00 to -20.00D

+10.00 to +2.50 in 0.50D steps

+2.00 to -8.00 in 0.25D steps

-8.50 to -20.00 in -0.50D steps

#### Vault Values

50, 100, 150, 200, 250, 300, 350, 400, 450, 500, 550

#### Wear & Replacement

Daily Wear. Replace at 6 months.

### Technical Support & Ordering

#### USA & Canada

Monday-Friday: 9:00AM to 8:00PM (EST)

#### Technical Consultation

Tel: 877.733.2012, option 2

E-Mail: [consultation@synergeyes.com](mailto:consultation@synergeyes.com)

#### Place Orders & Customer Care

Tel: 877.733.2012, option 1

Fax: 877.329.2012

E-Mail: [orders@synergeyes.com](mailto:orders@synergeyes.com)

Online: [www.synergeyes.com](http://www.synergeyes.com)

#### Outside USA & Canada

To locate an International Distributor visit [www.synergeyes.com/international](http://www.synergeyes.com/international).

To contact Customer Care in USA call +1 760.476.9410, option 1, or to place an order email: [intorders@synergeyes.com](mailto:intorders@synergeyes.com).