

KaVo ESTETICA® E70/E80 Vision

Get in touch with your vision.



KaVo. Dental Excellence.

Your wishes, our inspiration.

In the beginning there was an idea. A vision.

The development of a new treatment concept to provide optimum support in all work steps. A concept directed towards your requirements.

Firstly we listened intensively to what you had to say. With many of your suggestions in mind, we started on the path to implement exactly what you had desired: simpler operation, smooth work processes, modern patient communication, access to patient data, automated hygiene processes.

You requested us to retain tried and tested features – such as the unique ergonomic swivel chair concept.

It became our vision to develop our best treatment concept. A concept which fulfils the wishes of everyone – to make the daily routine easier.

The KaVo ESTETICA E70/E80 Vision – experience your vision.

Comfort according to individual requirements: The dentist element of the KaVo ESTETICA E70 und E80 Vision is naturally available again as a table, swing or cart version.*

** Available for KaVo ESTETICA E70 Vision*



All the innovations at a glance:

The new dentist element.

A simple touch.

- 1 Completely new operating concept for quick and direct access to all important functions for time-saving and smooth treatment processes

Patient communication.

Ready to go.

- 2 KaVo ERGOcam One, the easily operated intraoral camera
- 3 New high resolution screens: KaVo Screen HD in 22" and KaVo Screen One in 19"

CONEXIO.

Less searching. More finding.

- 4 System software CONEXIO, the smart data interface in your practice, supports you with useful functions and offers access to relevant patient data at all times

Hygiene.

Automated. Integrated.

- 5 Hygiene centre with integrated, removable adapters for instrument and suction hoses for time-saving, automatic cleaning and disinfection

New design.

Visibly modern. Noticeably gentle.

- 6 Blue Line Edition in ocean blue or smoke blue
- 7 Revised design with device components in modern dental white

Suspended chair.

Tried and tested.

- 8 New, easily adjustable armrests
- 9 Permitted patient weight up to 180 kg

Surgery.

Unrestricted. Integrated.

- 10 New surgical function with special small and convenient light surgical motor and integrated saline solution pump*

* See illustration on page 12

Just trust your intuition.

A simple touch

With the new sensitive touch display, your KaVo ESTETICA E70/E80 Vision is more intuitive than ever before!

With the touch display and its self-explanatory operating logic, you can access all important functions and information without any problems. Before, during and after treatment.

The completely new design of the user interface was specially developed to aid dental workflow. It is clearly designed, logically structured and ensures that you will find it easy to use at all times.


High-quality glass interfaces on the dentist and assistant ensure quick cleaning and improved hygienic features.

Working modes personalised for you.

Your treatment unit should be as individual as you are. You should therefore customise the equipment of your KaVo ESTETICA E70/E80 Vision towards your personal requirements. Specialisations are no problem at all due to the flexible instrumentation. Equip your dentist and assistant element exactly as you and your team need.

Many functions can be individually adapted to the practitioner, e.g. the name of the practitioner, the types of treatment or the direct key allocation on the start screen.

Further supplementary equipment and extensions are possible at any time thanks to future-oriented interfaces.

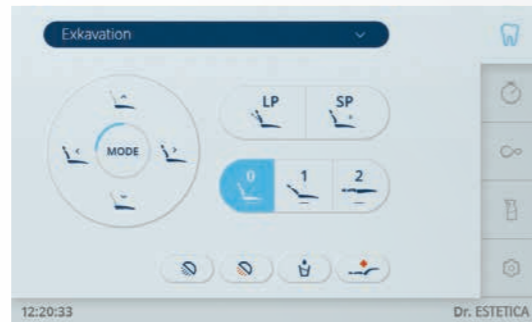

Well structured. Easy. Clear. The newly designed user interface of the dentist element.

Something special, not only for specialists.

More freedom. More comfort.

Conveniently access all functions and parameters of integrated special devices that you need during the treatment via touch display or foot control. You will maintain a constant view of all information on the dentist element. Individual treatment procedures can be freely programmed at any time.

Endodontics developed to a fine art.

Spare the acquisition of an expensive endodontic device and use the integrated Endo functions. You are well equipped for all endodontic applications with three different modes for arrow control.

You can set, call up and save all important parameters such as the transmission ratio and torque values.

The clever surgical option.

Make your work as comfortable as possible. You have the option to equip your KaVo ESTETICA E80 Vision with the new KaVo INTRA LUX S600 LED. Combined with the KaVo surgical instruments and the saline solution pump on the dentist element, the unrestricted functions of a modern surgical device will be at your disposal directly on your treatment unit.



The swing-arm version of the dentist element is also equipped with the new operating concept.



Due to the integrated and freely programmable Endo and surgical functions, your workflow will be complete, making the procurement of external devices unnecessary.



Everything for a healthy posture.

Make yourself comfortable.

Enjoy the overall ergonomic concept of the KaVo ESTETICA E70/E80 Vision. Comfortable sitting, standing and reclining – everything is possible.

There are practically no better prerequisites for a healthy posture and stress-free working for you and your assistants. Experience the freedom of natural, intuitive motion. And indulge your patients by providing comfortable reclining all the time. Comfortable upholstery and the height adjustment of the Comfort backrest ensures that your KaVo ESTETICA E70/E80 Vision can be adjusted to everybody's height.

Give your legs more freedom.

Enjoy the considerable legroom under the patient chair and position the foot control where you want it. The unique swing-chair concept of the KaVo ESTETICA E70/E80 Vision flexibility and comfort as you work. This means: a healthy posture, adequate positioning to the patient and pleasant freedom of movement in every treatment situation.



Maintain your posture!

Allow your natural movements to be supported and perform treatment with a healthy posture without physical stress. You can adjust the KaVo ESTETICA E70/E80 Vision to your personal needs and requirements. The variably adjustable height of the patient chair ranges from 350 mm to 900 mm. It guarantees optimal ergonomic conditions regardless of whether you sit or stand during treatment.

Treat your patients in the best position.

The patient positioning according to Trendelenburg is comfortable and ergonomic in every treatment position. Because the footrest of the patient chair is raised or lowered automatically to match the inclination of the backrest, the combined movement prevents the patient being squashed or stretched, while the armrest is being adjusted. The elevated leg position also supports the circulation system.

Indulge in some outstanding features.

The KaVo ESTETICA E70/E80 Vision and the PHYSIO Evo dentist chair, like their predecessors, were awarded the AGR (Campaign for Healthier Backs) seal of approval. Experts from different medical specialist fields have evaluated them, in accordance with stringent testing criteria, as an optimal ergonomic dental workplace.



PHYSIO Evo and Evo F offer perfect ergonomics with maximum comfort for healthy and fatigue-free working.



There is comfort in every detail.


Just right.

Offer your patients tailor-made comfort and be happy about the smooth flow of treatment in every situation. The KaVo ESTETICA E70/E80 Vision guarantees the full functional range of the treatment unit even for heavier patients, up to 180kgs. Neither you nor your patients need to forgo comfort.

Relaxation due to calm movement.

The slow motion function of your patient chair during treatment with the microscope means that it can be moved in slow motion. You will particularly appreciate this function when treating elderly patients.

*Lowest position 350 mm, highest position 900 mm.
Ideal for all requirements including working while standing.*

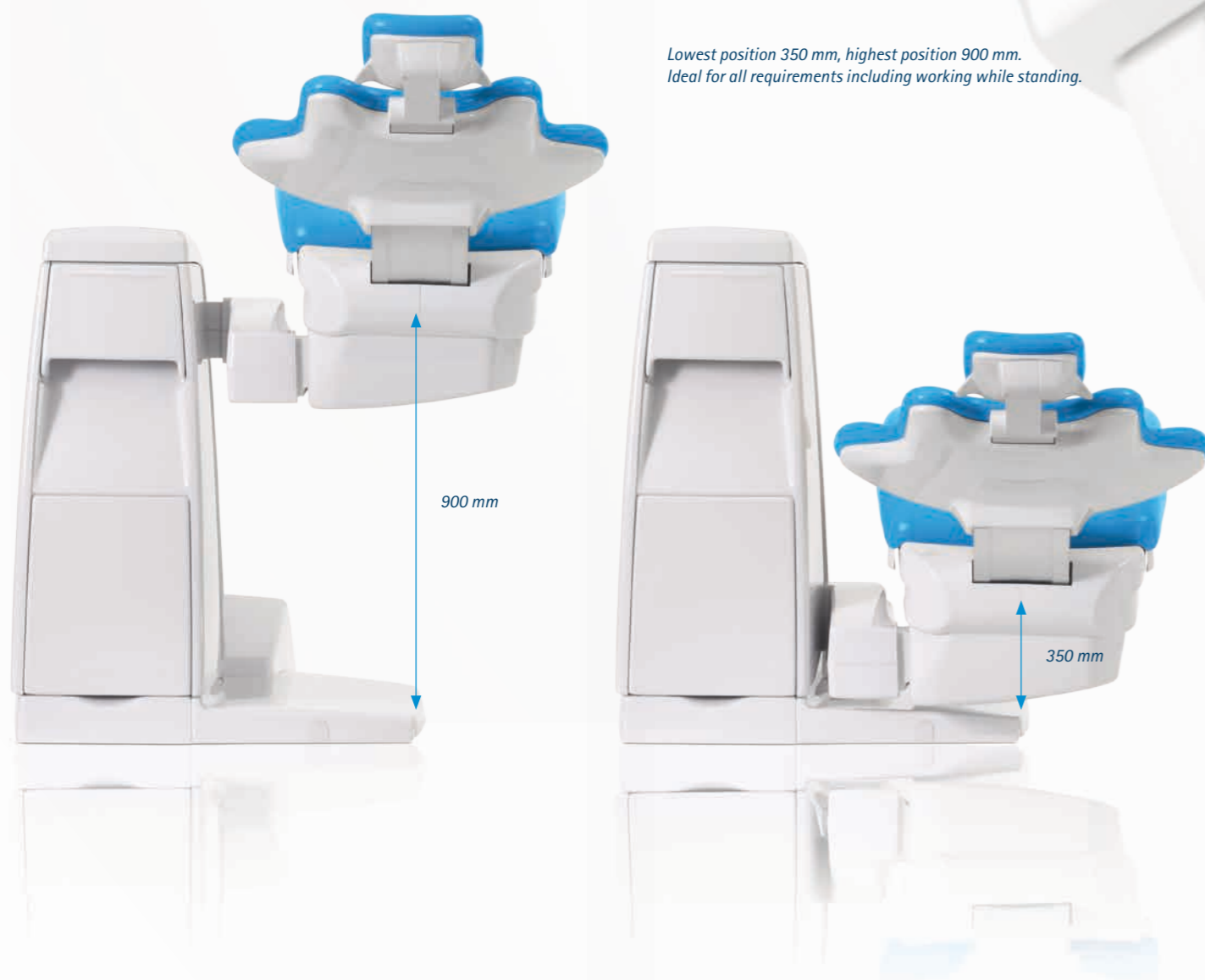

Easy to get into and out of the chair.

The new arm rests of the KaVo ESTETICA E70/E80 Vision make it easier for your patients to get into and out of the chair. They easily swing downwards and provide support and security during treatment, without you having to forgo freedom of movement.



The new arm rests not only look modern, they also provide additional support and security for you and your patient.

One starter for everything.

As your hands are occupied while treating your patients you can simply control your KaVo ESTETICA E70/E80 Vision with your foot. All functions and integrated devices can be used via the operating logic of the foot control. This also minimises the risk of cross contamination.

With the pleasant right-left motion of the foot control, your leg muscles will remain absolutely relaxed even during prolonged treatment. You can therefore work precisely and without fatigue. There is also a wireless foot control available as an option which provides additional flexibility in the positioning.



Perfect workflow means more efficient working patterns.

The perfect teamwork.

All instruments and functions are perfectly integrated in the KaVo ESTETICA E70/E80 Vision. Integrate endo-motors and surgical-motors without any problems. Control your dentist element directly via the touch display – or via the foot control. Expand your treatment unit per interface with USB-compatible devices and install these in a simple manner via the system software CONEXIO.

Enjoy maximum flexibility and superior performance with the dentist and assistant elements. The comfort is pioneering and the elements are quick, precise and efficient to use. Bringing ease and convenience to you and your team.

Everything finds its place.

Enjoy a perfect workflow. The KaVo Centro mounting system on which all important devices can be arranged sensibly and ergonomically in one place becomes the centre of your treatment processes, reaching paths and movement patterns. You have everything to hand and in easy reach.

You will also remain absolutely flexible if you wish to upgrade to further devices and future technologies. Devices and components are easily and quickly exchangeable due to the numerous adapters available.

The future will be easy – for the whole team.

The discovery of lightness.


Impress your team with the new assistant element of your KaVo ESTETICA E70/E80 Vision with its sleek and modern design. Intuitive and self-explanatory in operation. Equipped with a sensitive touch panel and the newly developed clear user interface. Comparable to the dentist element. All central chair and hygiene functions as well as the functions of operating lights are quickly and directly available.

Perfection for individual requirements.


Now it becomes personal. The assistant element is freely configurable and can be perfectly adapted to the actual requirements of your assistants. From the polymerisation handpiece to a supplementary saliva ejector.

A question of attitude.

Instruments can easily be precision adjusted to individual preferences for smooth and fast movement patterns. The removable handpiece tray provides more hygiene and the flexible height adjustment and the extensive swivel range of the assistant element provides pure ergonomics.

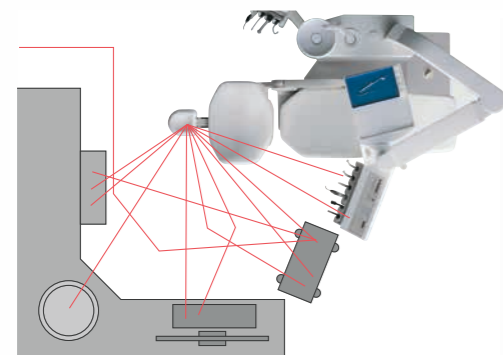

Movements with a conventional treatment centre.

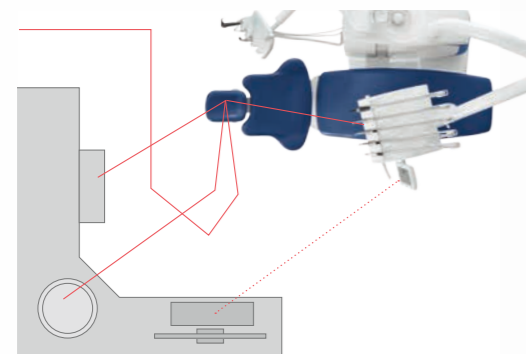

Movements with the new KaVo ESTETICA E80 Vision.

Customisable for individualised working – the Centro mounting system.



Precise positioning and stability – exclusively for the KaVo ESTETICA® E80 Vision.

Optimal workflow down to the smallest detail.

Experience a completely new level of working. Discover maximum treatment comfort for you and your patients. You will be delighted with the outstanding features available on the KaVo ESTETICA E80 Vision.

Stable and precise.

You are able to position your dentist element where you want to have it. The dentist element of your KaVo ESTETICA E80 Vision can be moved more precisely and positioned more stably with the serial pneumatic brake. So you will now have everything to hand – and everything under control.



Everything exactly where it must be.


Regardless of the height of your patient, the horizontal adjustment ensures that your working field – the mouth of the patient – will always stay in the same position. Your reaching distances will always be the same even if you move the patient chair while changing between upper jaw and lower jaw treatment.

The best position for every patient.

The chair can adapt perfectly to every patient. The motorised seat elevator means that patients with a hollow or round back and even children will be reclining comfortably and ergonomically. The seat elevator provides more reclining comfort, even during prolonged treatment, as you can reposition the patient as required.



You and your patients can enjoy the comfort of the motorised seat elevation system, exclusively offered by KaVo.

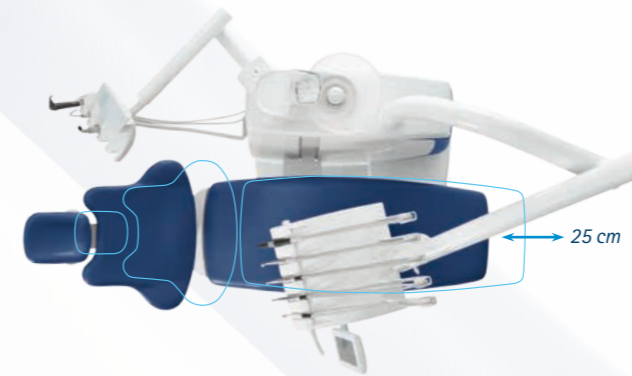

Use the horizontal movement as required, for example to provide space.

The spittoon bowl with automatic swing.

The spittoon bowl pivots automatically towards the patient and back, as soon the original treatment position is resumed – this offers added convenience.

Future proof – thanks to integrated USB interfaces.

Simply connect additional selected devices via the USB interface of your KaVo ESTETICA E80 Vision. This is really easy and your treatment system automatically recognises what has been connected to it.

NEW! Surgical motor S600 LED – powerful but light.

Work with a world champion! The new KaVo surgical motor S600 LED is one of the smallest and lightest surgical motors in the world. 28 % lighter and 20 % shorter than its predecessor model.

Surgery complete with integrated saline solution pump.

With the saline solution pump, directly integrated into the dentist element and the saline solution holder, you will experience maximum convenience and the most effective relief during every surgical treatment. It cannot get better.

With the KaVo S600 LED, one of the world's lightest and smallest micro motors, everyday surgical procedures become less strenuous.



Optimum patient communication. You have it in your hand.

Utilise the power of images.

Explain the necessary treatment measures to your patients in detail. Use the high resolution diagnostic examinations and images in brilliant colour reproduction for the patient communication. The CONEXIO patient communication system from KaVo has been especially developed for this purpose.

KaVo ERGOcam One – small camera, great results.

Easily explain the actual dental situation to your patients with images with high sharpness and outstanding colour reproduction. The depth of sharpness is so high you will not need to focus. The easy to operate KaVo ERGOcam One delivers perfect images at all times. This is not only convenient, it also saves valuable time. The image resolution can be activated by hand or with the foot control.

KaVo DIAGNOcam – X-ray vision without X-ray radiation.

Make yourself a picture or a video of the dental situation, on the dentist element or with the foot control. With the KaVo DIAGNOcam, you will detect carious lesions at an early stage. Light waves replace harmful radiation. Detect caries simply and be rewarded by a smile from your patient.

Dental microscope Leica M320 – sharp to the smallest detail.

Experience a new dimension of vision. With the Leica Microscope M320 that was specially developed for dentistry. Due to the combination of Full HD imaging technology and modern LED lighting, sensationally clear, sharp and brilliant images can be created in an outstanding sharpness of depth. With the MultiFoc Objektiv and the variable focal distance of 200 to 300 mm, you can change between different enlargement levels and will always be in the right focus.

The microscope images are displayed on the KaVo Screen and can be conveniently saved in relation to the patient in CONEXIO by simply pressing a button.

You will see considerably more with the Leica M320 dental microscope.

Show it to your patients.

Impress your patients with sharp images. The new screens KaVo Screen One and KaVo Screen HD provide high-resolution image quality, convincing contrast values and brilliant colour reproduction. Optimally coordinated to all functions of imaging devices like the camera or microscope. The screens are easily operated via the touch display on the dentist element.

For maximum hygiene safety, both screens have glass safety panels, a narrow gap on the housing and can be adjusted with one hand.

KaVo Screen One – the introduction into the HD world

- HD-true flat screen
- aspect ratio: 16:9
- screen diagonal: 19 inch

KaVo Screen HD – the perfect choice in full HD

- Full HD-true flat screen
- aspect ratio: 16:9
- screen diagonal: 22 inch
- Two digital inputs for microscope and PC



CONEXIO joins what belongs together.

Whatever you intend to do, CONEXIO opens the way.

Experience more efficiency in your everyday routine with CONEXIO. Access all patient data directly via the touch display of the treatment unit. Need to add further data? No problem. All collected data are directly assigned to the patient file. CONEXIO makes workflows easier. Let yourself be guided logically and highly efficiently through all work processes. Time-saving. User-friendly. Easy.



Helpful for your diagnostics. Educational for your patients.

Use the practical CONEXIO screenshot option to cut out an image from a certain area and save as a separate image. Whether from a video or from an image – you will no longer need any further programs or alternative computers for trimming, contrast, brightness, enlargement and storage.

Maintain an overview – with the automatic diagnostic history.

Every diagnosis will be well logged, summarised and saved, completely automatically. Use search filters to find specific information quickly in the history.

A treatment with reliability.

Your work will be a pleasure. Everything is exactly where you want and need it, namely on your dentist element. The user-friendly interface on the touch display of the treatment unit makes this possible and guarantees full access to the patient file. And at times when you have no free hand, all central CONEXIO functions can also be controlled and activated via the foot control. Convenient and hygienic.

Always up-to-date.

You would like to connect a new treatment device to your KaVo ESTETICA E70/E80 Vision? There is nothing easier! CONEXIO recognises KaVo hardware components immediately and automatically. This means: Plug and Play for every system expansion.


Soft and hardware updates are recognised automatically and are prepared for installation. Thus your system is always up to date and the security of your data inventory assured.

Swap Tray – Your protected exchange platform.

With Swap Tray, data such as digital photos, X-ray pictures or videos can be collected and shared with other CONEXIO users. Images can be exchanged in standard format and also anonymised on request. It's that easy and without burdening the email box.

Secure is secure – with automatic data backup.

Safe can be so easy. Especially if you don't have to think of the backup any more. You decide when and exact data backups are generated by CONEXIO periodically. Thus you will have more time for the essential matters.

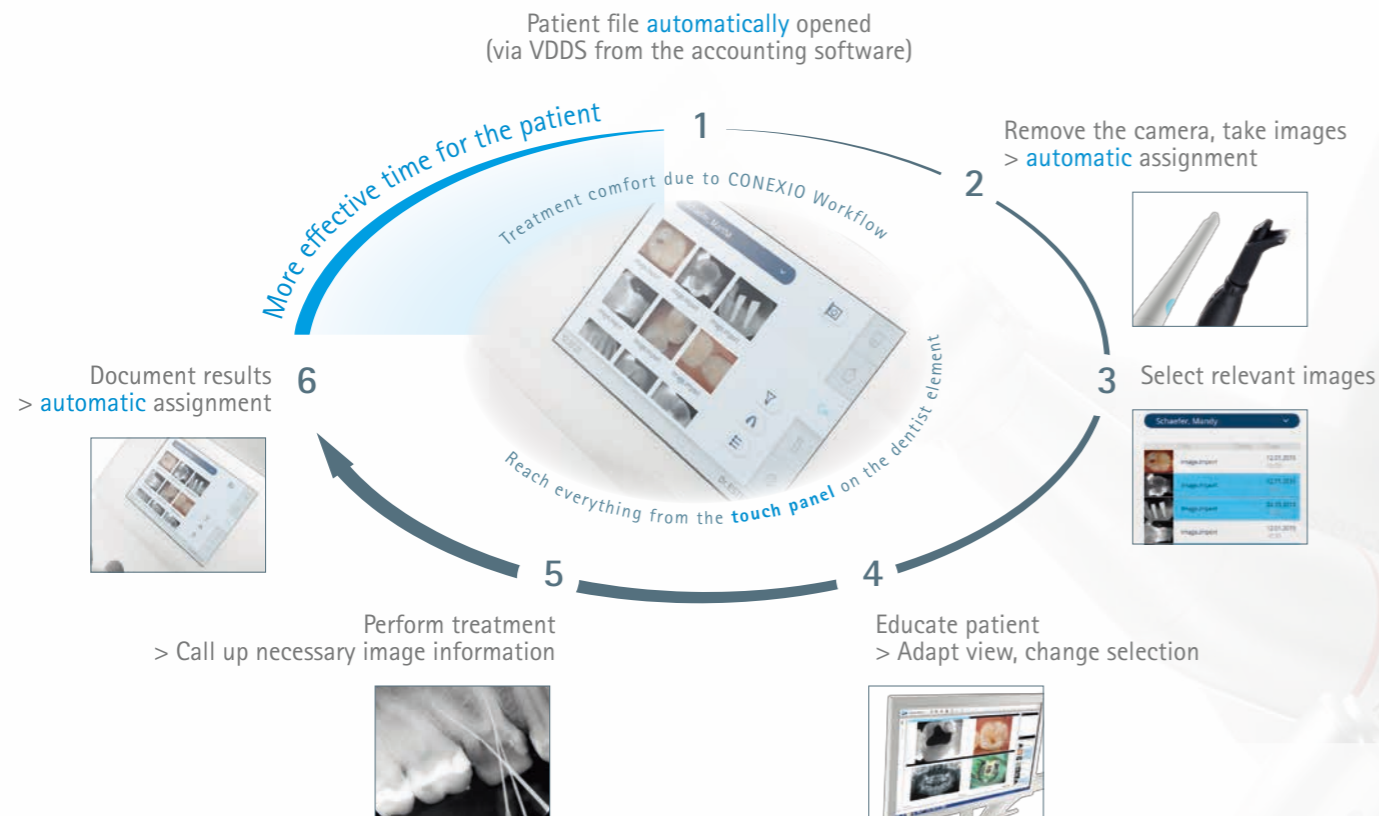

Fascination in form and colour.

Our vision for style and design.

View the KaVo ESTETICA E70/E80 Vision from a completely different perspective. The colour and form focus on style and class. It is no wonder that the KaVo ESTETICA E80, designed by SLOGDESIGN, was awarded the renowned Good Design Award. Elegant unit body. High-quality surfaces. Ergonomic overall design. And numerous opportunities for individualisation due to the great selection of paint and upholstery colours. The entire appearance invokes trust and gives a feeling of quality. A vision that you, your patients and your entire team feel good about.

Metallic paint colours



Plain colours



Upholstery colours



Comfort backrest: maximum comfort and wide shoulder rest for patients.

Perfect position for your patients.

All patients are different so it is essential to have a patient chair which takes this into consideration.

Your patients will recline very comfortably with the wide Comfort backrest. Thanks to its height adjustment system, the KaVo ESTETICA E70/E80 Vision adapts to every body height and constitution. The Progress backrest is especially narrow for optimum access to patients. This special head cushion allows your patients to recline comfortably even during long treatments. Ergonomics at the highest level and comfort too.

Comfort begins at the head.

Make yourself and your patients comfortable. Adjust the headrest of your KaVo ESTETICA E70/E80 Vision exactly to the position of the head and mouth. The inclination angle of the head is infinitely adjustable. The fixed positions are precisely adapted to the anatomy. Whether motorised, pneumatic or mechanical, you will always find the optimum position for the headrests. Comfortable and precise.

Progress backrest: correct working posture and optimum access for the dentist, especially when treating children.

The double-jointed headrest is available in the pneumatic, motorised and mechanical version.



When practical intelligence turns into safe hygiene.

Ready for operation. With safety.

Do not let your time be wasted with a complicated hygiene routine. You can be sure, the integrated hygiene functions and the fully automatic processes of your KaVo ESTETICA E70/E80 Vision will look after this, highly efficiently and absolutely effectively. You will be guided in the operation of the hygiene functions via the touch display. A real time saver.

Automated, integrated.

Rely on something new with the simple and reliable automatic cleaning and disinfection of the suction and instrument hoses. The hygiene centre is now newly integrated in the unit body of the KaVo ESTETICA E70/E80 Vision. You can connect both cleaning and instrument hoses here. The hygiene centre can be removed and thermodisinfected.

Pleasant. Removable parts.

Typically, KaVo has also given great thought to the areas that your hands touch. You can simply remove contaminated device parts such as handles, instrument trays or spittoon bowls. They can be cleaned, disinfected and some even sterilised, depending upon circumstances.



Handles, silicone mats, spittoon bowls, instrument trays and swing-arm clamps may be easily removed and disinfected.

Improved protection with an anti-reflux system.

Protect your water supply system with the anti-reflux technology in KaVo instruments and motors. It prevents contaminated spray water from being aspirated into the treatment centre.

Automatic rinsing program.

Simplify the daily rinsing of all water conveying lines with the automatic cleaning functions. Four programs guarantee reliable hygiene at the press of a button.

During normal operation, continuous germ reduction with OXYGEN-AL 6, prevents microorganisms from multiplying and spreading. The impeccable quality of the cooling and rinsing fluids guarantees an intensive germ reduction ("weekly" cleaning program) even during extended downtimes.

HYDRoclean cleans and disinfects the entire suction and discharge system and the suction hoses. Furthermore, DEKASEPTOL Gel is dosed automatically, eliminating the manual process for assistants.



OXYmat und DEKAmat, both integrated into the device - now introduce standardised and automated hygiene into your practice. Furthermore, automated care promotes operational reliability.



Experience efficiency with the right equipment.

Ready for the future, today.

Work in the future, today. Because with the wide equipment range of the KaVo ESTETICA E70/E80 Vision, you are optimally prepared for everything that comes as well as everything that goes.

Experience your requirements.

Only accept the best in instrumentation. The new KaVo MASTER-matic series covers all application areas and you will work with maximum quality and precision. Enjoy supreme vision and optimum access. The M25L quick contra-angle handpiece with a small head size offers improved access and visibility.

Be prepared for the future.

With the integrated USB interfaces of the KaVo ESTETICA E70/E80 Vision, your treatment unit can be extended by further components at any time. In this way, you will have more freedom in the future with the procurement of supplementary devices.

Put your practice in the best light.

Rely on the KaVoLUX 540 LED 5-star light. A natural white light of highest quality. The light can be switched on and off and operated without contact, even in COM-POsave mode. Closed housing with smooth interfaces and removable, thermo-disinfectable handles provide quick and reliable hygiene.

Optimum balance of power and output.

PiezoLED, the high-precision device for thorough teeth cleaning.

With the PiezoLED Ultrasonic scaler, professional tooth cleaning is not only precise, it also saves your valuable time. The titanium sonotrodes of the PiezoLED ensure controlled linear oscillations, enabling efficient treatment of your patients. Due to the precise operation, teeth and gingiva remain undamaged even in delicate areas.

The circular light emissions produce optimum viewing conditions, and in combination with the lightweight and ergonomic design you can always treat your patients under optimum conditions.


Micromotor INTRA LUX KL 703 LED – something so great can be small.

Experience what a small motor can achieve. Endless power. Fatigue-free working. Speed range between 100 and 40.000 min⁻¹ with a 1:1 transmission ratio. Easy. Compact design. Pleasant in your hand. SMARTdrive-optimised, you can work highly precisely and efficiently with this micromotor in every treatment situation.

Convince yourself of its outstandingly smooth running which makes precise working much easier and convenient.

Full power: SMARTdrive.

Power inside! Rely on precise, efficient and time-saving work. With the SMARTdrive motor control. It displays its torque even in low speed ranges of 100 min⁻¹. Frequent refitting of the contra-angle handpiece for the removal of caries and polishing is unnecessary. You can perform the treatment with only one contra-angle handpiece and not waste any valuable time. You are in the pole position in every speed range. You need only accelerate.

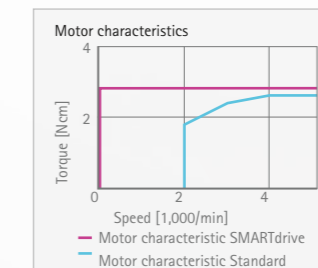

Leverage the high torque and vibration-free start-up characteristics of your KaVo SMARTdrive with its extended speed range of 100 - 40,000 rpm.

PiezoLED: the innovative ultrasonic technology for the highest standard of patient treatment.

Equipment options.

	E80 Vision		E70 Vision				E80 Vision		E70 Vision		
	T	C	T	S	C		T	C	T	S	C
Patient chair											
Soft upholstery	●	●	●	●	●		●	●	●	●	●
Mot. Horizontal tracking	●	●	-	-	-		●	●	●	●	●
Mot. seat elevation	●	●	-	-	-		●	●	●	●	●
Backrest											
Comfort	x	x	x	x	x		○	○	○	○	○
Progress	x	x	x	x	x		○	○	○	○	○
Headrest											
Double jointed mechanical	-	-	●	●	●		○	○	○	○	○
Double-jointed pneumatic	●	●	○	○	○		○	○	○	○	○
SOMATIC E-motorised	○	○	○	○	○		○	○	○	○	○
Comfort cushion	○	○	○	○	○		○	○	○	○	○
Armrest											
Left	○	○	○	○	○		●	●	●	●	●
Left and right	○	○	○	○	○		●	●	●	●	●
Dentist element											
Touch display	●	●	●	●	●		●	●	●	●	●
1. Turbine hose	●	●	●	●	●		●	●	●	●	●
2.-5. Turbine hose	○	○	○	○	○		○	○	○	○	○
1. KL motor 703	●	●	●	●	●		●	●	●	●	●
2.-5. KL motor 703	○	○	○	○	○		○	○	○	○	○
Endo function	○	○	○	○	○		○	○	○	○	○
COMFORTbase	○	○	○	○	○		○	○	○	○	○
Triple function handpiece	-	-	●	●	●		●	●	●	●	●
Multifunctional handpiece	●	●	○	○	○		●	●	●	●	●
PiezoLED ultrasonic scaler	○	○	○	○	○		○	○	○	○	○
Physiological saline solution	○	○	○	○	○		○	○	○	○	○
Surgical motor S600 LED	○	○	-	○	-		●	●	●	●	●
Spray heating for instruments	●	●	●	●	●		○	○	○	○	○
Anti-reflux valve for instruments	●	●	●	●	●		○	○	○	○	○
Removable handpiece tray	●	●	●	●	●		○	○	○	○	○
Small X-ray viewer	○	○	○	○	○		○	○	○	○	○
Large X-ray viewer	○	○	○	○	○		○	○	○	○	○
USB interface	●	●	○	○	○		○	○	○	○	○
Pneumatic brake for arm system	●	●	○	○	○		○	○	○	○	○
Tray holder for 2 standard trays	●	●	●	●	●		○	○	○	○	○
Adjustable service table for two standard trays	●	-	●	-	-		○	○	○	○	○
Service table 1658	-	●	-	-	●		○	○	○	○	○
Assistant element											
Glass touch panel							●	●	●	●	●
Spray mist suction							●	●	●	●	●
Saliva ejector							●	●	●	●	●
2nd Saliva ejector							○	○	○	○	○
Triple function handpiece							○	○	○	○	○
Polymerisation handpiece Satelec Mini LED							○	○	○	○	○
Support for tray holder							○	○	○	○	○
Height-adjustable arm							●	●	●	●	●
Removable handpiece tray							●	●	●	●	●
Unit Body											
Glass spittoon bowl							x	x	x	x	x
Porcelain spittoon bowl							x	x	x	x	x
Motor-driven, pivoting spittoon							●	●	○	○	○
Hygiene centre							●	●	●	●	●
Permanent germ reduction							●	●	●	●	●
Intensive disinfection							●	●	●	●	●
Hydroclean (rinsing system)							●	●	●	●	●
DEKAmat and OXYmat							●	●	●	●	●
Warm water boiler							●	●	●	●	●
VACUstopp							●	●	●	●	●
Wireless foot control							○	○	○	○	○
Supply system											
Water block DVGW							●	●	●	●	●
Additional device connections							○	○	○	○	○
AQUAmat							●	●	●	●	●
CENTRAmat							●	●	●	●	●
Disposal system											
External suction							●	●	●	●	●
DÜRR amalgam separator							○	○	○	○	○
Dürr automatic separation unit							○	○	○	○	○
Operating light											
KaVoLUX 540 LED							○	○	○	○	○
Operating light mounting post							○	○	○	○	○
Patient communication											
KaVo Screen One							○	○	○	○	○
KaVo Screen HD							○	○	○	○	○
CONEXIO							○	○	○	○	○
KaVo ERGOcam One intraoral camera							○	○	○	○	○
KaVo DIAGNOcam							○	○	○	○	○

● Serial equipment ○ Optional equipment x Must be selected - Not available



KaVo. Dental Excellence.