



**clearlab**®

## GIVE YOUR PATIENTS THE BEST WITH CLEAR 55 A

Optimized Comfort and Vision for all Myopes and Hyperopes

### Improved Comfort:

- Extended Wearing Hours with AQUAGRIP Technology
- Superior Comfort with TrueForm Edge Technology
- Excellent Fit with superior design and geometry

### Improved Vision:

- A high-definition vision with aspheric and advanced optics
- Enhanced Depth of Focus & better night vision with larger optic zone
- A healthy & active lifestyle with UV protection in the lens

### Improved Packaging:

- Small and neater packaging
- Easy to open
- Very convenient to store and carry

Optimizing Visual Acuity, Wearer Comfort, Fit and Handling:



The peripheral and edge design of clear55A provides a balance between excellent comfort and assured handling.

The larger diameter of clear55A assists in centering the lens on the eye and ensures a wide latitude of fit over a wide range of corneal sizes, providing practitioners with a versatile fitting lens



The visual acuity of clear55A is optimized over the entire power range and the incorporation of aberration correcting aspherics provide wearers with maximum contrast sensitivity, particularly in low illumination where it is needed the most.



clear55A provides an optimal balance of handling, comfort and vision.

The lens has been enhanced with aberration correcting aspheric optics to optimize the overall visual experience and encompasses a design philosophy that addresses the wearer's handling and comfort needs.

For healthy active lifestyle clear55A provides UV protection.




# clear 55A

The Clear55A range from Clearlab® offers its users an excellent high-definition vision with an advanced optics, along with a superior long-lasting comfort and excellent handling by keeping a stable moisture circulation throughout a long wearing schedule for a better health of the cornea, and thus was achieved through different innovations in the manufacturing process which are:


## AQUAGRIP technology

It keeps the lens hydrated all day long reducing dropout ratios and promoting comfort throughout the wearing hours.

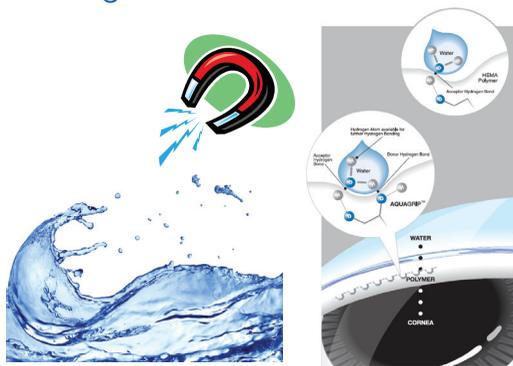


## Truform EdgeTechnology

Promotes ultra-thin shape with rounded edge design to keep proper tear film exchange and overall easy lens handling.

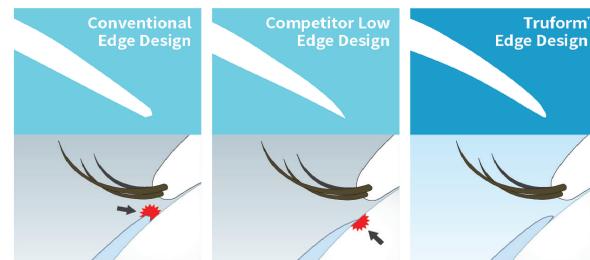

## High-Definition Multi-Band Aspheric

Enhanced depth of focus due to the enlarged optic zone with an advanced aspheric optics.

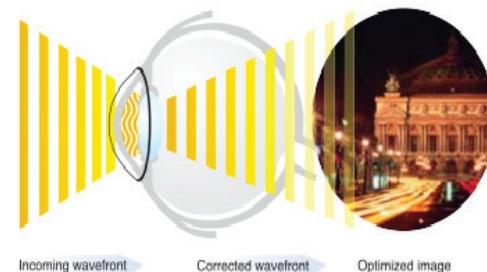

## Product Specifications:

Name	Clear55A™ (Soft, Clear and HD Aspheric Contact Lenses)
Material	Hydrogel lens with 45% Methafilcon A & 55% water content
Wearing Modality	Monthly Lens with daily wear
Base Curve	8.70mm
Diameter	14.50mm
Oxygen Permeability	23.4 DK
Center Thickness	0.09 mm @ -3.00 D
Handling Tint	Light Blue
Oxygen Permeability DK@35oC	23.43 X 10 <sup>-11</sup> (cm <sup>2</sup> /s) (mlO <sub>2</sub> /ml.mm Hg)
UV Blocking	UVA + UVB
Power Range	-0.25 to -5.00 (0.25 steps) -5.50 to -10.00 (0.50 steps) +0.50 to +6.00 (0.50 steps)
Package size	6 lenses per pack